

# Quick Start Guide

LinkSquare Basic

---


This Quick Start Guide contains the **UDID**  
for your LinkSquare Basic Software.

Learn more at [linksquare.io/help](https://linksquare.io/help)



LinkSquare

# Product Diagram

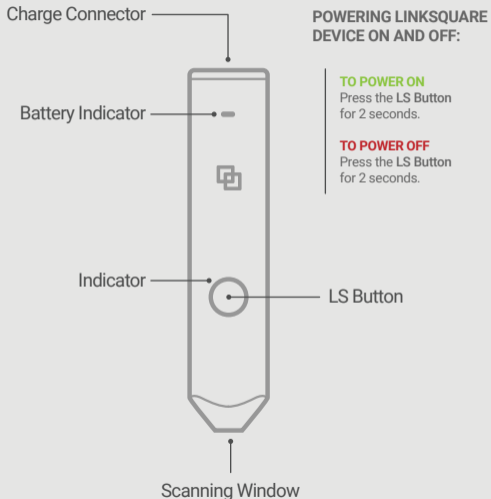

# More About LinkSquare



Before using your LinkSquare, please read carefully all the safety precautions in the **Product and Safety Instructions**.



## REGULATORY INFORMATION

See the list of country certifications in the **Product and Safety Instructions**.



## CUSTOMER SUPPORT


Need help?  
Want more information?  
[linksquare@stratiotechnology.com](mailto:linksquare@stratiotechnology.com).



## JOIN THE LINKSQUARE MOVEMENT

 [fb.com/LinkSquareIR](https://fb.com/LinkSquareIR)  
 [twitter.com/LinkSquareIR](https://twitter.com/LinkSquareIR)


# Getting Started

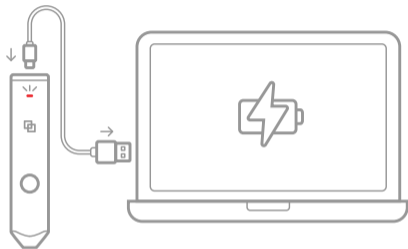

1

## FULLY CHARGE YOUR LINKSQUARE DEVICE PRIOR TO USE.

The **Battery Indicator** will stop blinking when charging is complete.

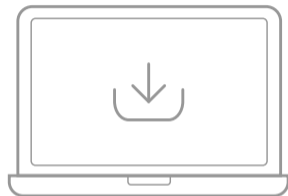


2

## DOWNLOAD LINKSQUARE BASIC ON YOUR COMPUTER.

Download the software from the link provided in your order-related emails.

\*Windows 7 or higher is required to install.

# Getting Started

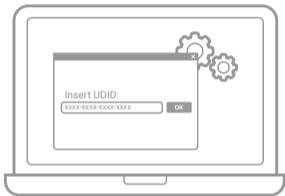

3

**ENTER YOUR UDID AND  
INSTALL THE SOFTWARE.**

Enter the following UDID:



4

**TURN ON YOUR  
LINKSQUARE DEVICE.**

The **Indicator** turns blue when LinkSquare is on.

# Connecting to Computer

There are two options to connect LinkSquare device and your computer.  
Press the **LS Button** to switch between the two Modes.

## IoT MODE

### CONNECT VIA YOUR INTERNET ROUTER.

Ensure that the **Indicator** is set to **blue** for IoT Mode.



Select "**IoT Mode**" icon on the launch page of the software and follow the directions.

## AP MODE

### CONNECT VIA LINKSQUARE'S BUILT-IN WI-FI.

Ensure that the **Indicator** is set to **green** for AP Mode.



Search for "LS-XXXXXXX" Wi-Fi network on your computer and enter the following password: **00000000**.

# Playing with LinkSquare Basic

1




Connect LinkSquare device and LinkSquare Basic software by pressing "Connect" icon in the software. The **Indicator** turns white when LinkSquare is connected to the software.

2



Point LinkSquare device straight down at the object to be scanned and press **LS Button** or click "Scan" icon in the software to scan.

3

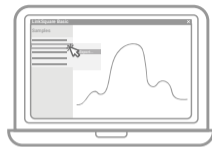

You can view the collected data by selecting one from the "Sample" window. Right click on the sample for more options. You can also export the collected data to various output file formats.

4



Click "Disconnect" icon in the software to disconnect LinkSquare from the software.

# Exploring LinkSquare App



- Check out our mobile app made from
- LinkSquare Basic. You can explore the
- included Color Palette app and see when
- new apps are released.

• **Note that LinkSquare Basic does not include**

• **a mobile app builder.**

## 1 DOWNLOAD AND LAUNCH THE LINKSQUARE APP.



Currently supporting  
iOS versions only.  
**iOS 10 or higher is  
required to install**



Android version  
is coming soon!

## 2 TURN ON LINKSQUARE DEVICE AND CONNECT TO YOUR MOBILE DEVICE.

Please refer to **Connecting to Computer** section for two Modes of connection.

## 3 CONNECT LINKSQUARE DEVICE AND LINKSQUARE APP.

Press "Connect" on the app home screen and the **Indicator** will turn white when connected.

## 4 HAVE FUN EXPLORING THE LINKSQUARE APP!